

**产品名称: UDG Rabbit Polyclonal Antibody**

**产品货号: APRab19603**

### 产品概述 (Summary)

描述(Description)	Rabbit polyclonal Antibody
宿主(Host)	Rabbit
应用(Application)	WB,IHC,ICC/IF,ELISA
种属反应(Reactivity)	Human,Mouse,Rat
偶联物(Conjugation)	Unconjugated
修饰(Modification)	Unmodified
同种型(Isotype)	IgG
克隆(Clonality)	Polyclonal
剂型(Form)	Liquid
保存条件(Storage)	Store at 4°C short term. Aliquot and store at -20°C long term. Avoid freeze/thaw cycles.
储存溶液(Buffer)	Liquid in PBS containing 50% glycerol, 0.5% protective protein and 0.02% New type preservative N.
纯化 (Purification)	Affinity purification

### 产品应用(Application)

稀释比(Dilution Ratio)	WB 1:500-1:2000,IHC 1:100-1:300,ICC/IF 1:50-1:200,ELISA 1:20000-1:40000
分子量(Molecular Weight)	35kDa

### 抗原信息(Antigen Information)

基因名(Gene Name)	UNG
别名(Alternative Names)	UNG; DGU; UNG1; UNG15; Uracil-DNA glycosylase; UDG
基因 ID(Gene ID)	7374.0
SwissProt ID	P13051
免疫原(Immunogen)	The antiserum was produced against synthesized peptide derived from human UNG. AA range:191-240

### 研究背景 (Background)

This gene encodes one of several uracil-DNA glycosylases. One important function of uracil-DNA glycosylases is to prevent mutagenesis by eliminating uracil from DNA molecules by cleaving the N-glycosylic bond and initiating the base-excision

repair (BER) pathway. Uracil bases occur from cytosine deamination or misincorporation of dUMP residues. Alternative promoter usage and splicing of this gene leads to two different isoforms: the mitochondrial UNG1 and the nuclear UNG2. The UNG2 term was used as a previous symbol for the CCNO gene (GeneID 10309), which has been confused with this gene, in the literature and some databases. [provided by RefSeq, Nov 2010],disease:Defects in UNG are a cause of immunodeficiency with hyper-IgM type 5 syndrome (HIGM5) [MIM:608106]. Hyper-IgM syndrome is a condition characterized by normal or increased serum IgM concentrations associated with low or absent serum IgG, IgA, and IgE concentrations. HIGM5 is associated with profound impairment in immunoglobulin (Ig) class-switch recombination (CSR) at a DNA precleavage step.,function:Excises uracil residues from the DNA which can arise as a result of misincorporation of dUMP residues by DNA polymerase or due to deamination of cytosine.,online information:UNG mutation db,PTM:Isoform 1 is processed by cleavage of a transit peptide.,similarity:Belongs to the uracil-DNA glycosylase family.,subunit:Monomer. Interacts with HIV-1 Vpr.,tissue specificity:Isoform 1 is widely expressed with the highest expression in skeletal muscle, heart and testicles. Isoform 2 has the highest expression levels in tissues containing proliferating cells.,


### 研究领域 (Research Area)

Base excision repair;Primary immunodeficiency;


### 注意事项 (Note)

For research use only.


### 图片 (Image Data)

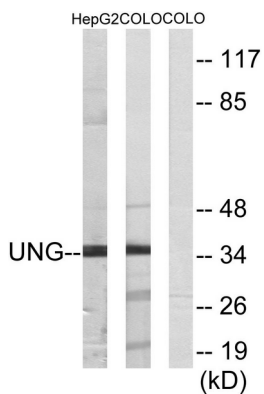


Western blot analysis of lysates from HepG2 and COLO cells, using UNG Antibody. The lane on the right is blocked with the synthesized peptide.

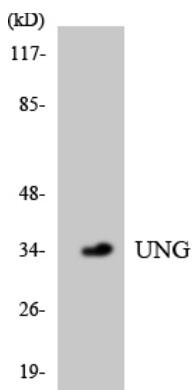

Western blot analysis of the lysates from HepG2 cells using UNG antibody.

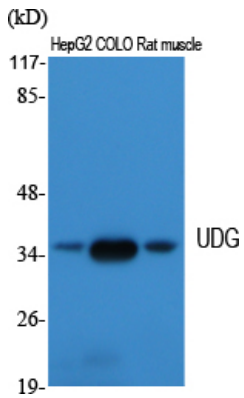

Western Blot analysis of various cells using UDG Polyclonal Antibody. Secondary antibody was diluted at 1:20000

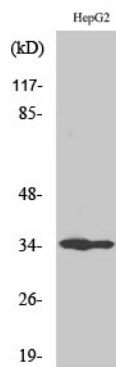

Western Blot analysis of COLO205 cells using UDG Polyclonal Antibody. Secondary antibody was diluted at 1:20000

## EnkiLife 优势产品推荐

产品类别	产品货号	产品名称
WB 解决方案	RA10020	2h 极速 WB 即用型全流程试剂盒
	RA10021	4h 快速 WB 即用型全流程试剂盒
	RA10042	5h 畅享版 WB 全流程试剂盒
	RA10037	校准级彩色预染蛋白 Marker (8-180kDa)
	RA10038	校准级彩色预染蛋白 Marker (10-250kDa)
	RA10039	校准级高分子彩色预染蛋白 Marker (25-400kDa)
TSA 多重荧光染色试剂盒	RA10008	TSA 双标三色多重荧光染色试剂盒 (mIHC)
	RA10009	TSA 三标四色多重荧光染色试剂盒 (mIHC)
	RA10010	TSA 四标五色多重荧光染色试剂盒 (mIHC)
	RA10011	TSA 五标六色多重荧光染色试剂盒 (mIHC)
	RA10012	TSA 六标七色多重荧光染色试剂盒 (mIHC)
IHC 检测试剂盒	RA10006	HRP Anti-Mouse/Rabbit IHC Detection System
	RA10007	Polymer-HRP Anti-Mouse/Rabbit IHC Detection System
抗体标记试剂盒	RE80004p	辣根过氧化物酶(HRP)抗体标记试剂盒
	RE80002q	Sulfo-NHS-生物素标记试剂盒
	RE80007p	Cy3 荧光素标记试剂盒
	RE80011p	Fluor488 荧光素标记试剂盒
	RE80017p	Fluor750 荧光素标记试剂盒
	RE80005p	藻红蛋白(R-PE) 抗体快速标记试剂盒
	RE80040	PE-Cy7 串联染料抗体快速标记试剂盒
稳转细胞系构建服务 (免费赠送全膜 WB 验证)	TS-0001	过表达稳转细胞系构建
	TS-0002	敲低稳转细胞系构建
	TS-0003	敲除细胞系构建

限于篇幅，更多产品详见官网：<https://www.enkilife.cn>