


概述 (Summary)

英文全称	TIMP-2/Metalloproteinase inhibitor 2/Tissue Inhibitor of Metalloproteinases 2
纯度 (Purity)	Greater than 95% as determined by reducing SDS-PAGE
内毒素 (Endotoxin level)	<1 EU/μg as determined by LAL test.
蛋白构建 (Construction)	Recombinant Human Tissue Inhibitor of Metalloproteinases 2 is produced by our Mammalian expression system and the target gene encoding Cys27-Pro220 is expressed with a 6His tag at the C-terminus.
Accession #	P16035
蛋白标签 (Tag)	
表达宿主 (Host)	Human Cells
种属 (Species)	Human
预测分子量 (Predicted MW)	22.79 KDa
蛋白形态 (Form)	Lyophilized from a 0.2 μm filtered solution of 20mM PB, 150mM NaCl, pH 7.2.
储存缓冲液 (Buffer)	
运输方式 (Shipping)	The product is shipped at ambient temperature. Upon receipt, store it immediately at the temperature listed below.
稳定性&储存 (Stability &Storage)	Store at ≤-70°C, stable for 6 months after receipt. Store at ≤-70°C, stable for 3 months under sterile conditions after opening. Please minimize freeze-thaw cycles.
复溶 (Reconstitution)	Always centrifuge tubes before opening. Do not mix by vortex or pipetting. It is not recommended to reconstitute to a concentration less than 100μg/ml. Dissolve the lyophilized protein in distilled water. Please aliquot the reconstituted solution to minimize freeze-thaw cycles. Always centrifuge tubes before opening. Do not mix by vortex or pipetting. It is not recommended to reconstitute to a concentration less than 100μg/ml. Dissolve the lyophilized protein in distilled water. Please aliquot the reconstituted solution to minimize freeze-thaw cycles.

电泳图 (SDS-PAGE image)

Product Name: Recombinant Human TIMP-2 (C-6His)
Catalog #: PHH1655

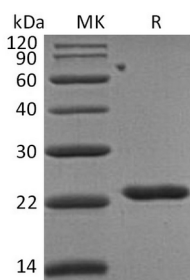

背景 (Background)

分子别名 (Alternative Names)

Metalloproteinase Inhibitor 2; CSC-21K; Tissue Inhibitor of Metalloproteinases 2; TIMP-2; TIMP2

背景介绍 (References)

Tissue inhibitors of metalloproteinases or TIMPs are a family of proteins that regulate the activation and proteolytic activity of the zinc enzymes known as matrix metalloproteinases (MMPs). There are four members of the family, TIMP-1, TIMP-2, TIMP-3, and TIMP-4. Tissue Inhibitor of Metalloproteinases 2 (TIMP-2) is a non N-glycosylated protein with a molecular mass of 22 kDa. It produced by a wide range of cell types, which inhibits MMPs non-covalently by the formation of binary complexes and irreversibly inactivates them by binding to their catalytic zinc cofactor. TIMP-2 also has erythroid-potentiating and cell growth promoting activities.

注意事项 (Note)

For Research Use Only , Not for Diagnostic Use.